



# Cellular ATMP Fresh Allogeneic Example Receipt Checklist FINAL

**Creator:** University Hospital of Wales & University Hospitals Bristol & Weston NHS Foundation Trust

**Document version number:** V1.0

**Date written:** November 2020


## End user rights:

This document is shared with permission for re-use to distribute, remix, adapt, and build upon the material in any medium or format for non-commercial purposes only, so long as the attributions listed below are given.

**Attributions:** University Hospital of Wales & University Hospitals Bristol & Weston NHS Foundation Trust

This document is made available under a Creative Commons Attribution - NonCommercial 4.0 International License as described here:  
<https://creativecommons.org/licenses/by-nc/4.0/>

The information, materials and any opinions contained in this document are provided for general information and educational purposes only, are not intended to constitute legal or other professional advice and should not be relied on or treated as a substitute for specific advice relevant to particular circumstances. Although we make all reasonable efforts to ensure the information is up to date, we make no representations, warranties or guarantees in that regard. In no event shall the creator(s) be liable for any direct, indirect, special, consequential or other damages that are related to the use or reliance whatsoever in the content of the document or any part thereof.

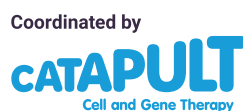

## Example checklist for receipt and storage of allogeneic, fresh, cellular ATMPs


**18 Nov 2020**

This document sets out an example generic checklist for receipt and storage of fresh, allogeneic, cellular Advanced Therapy Medicinal Products (ATMPs). It may require adaptation to accommodate local procedures or product-specific requirements. The checklist should be adjusted to reflect whether the ATMP is received direct from the manufacturer or via a Stem Cell Laboratory. Institutional readiness checklists produced by the Pan-UK Pharmacy ATMP Working Group for relevant cellular ATMP type should also be considered. (Available from <https://www.sps.nhs.uk/>).

Separate example checklists are also available for autologous fresh cellular ATMPs and cryopreserved autologous/allogeneic cellular ATMPs.

## FRESH ALLOGENEIC CELLULAR ATMP PRODUCT RECEIPT & STORAGE CHECKLIST

Product Name				
Sending laboratory if applicable				
Manufacturer (if different to above)				
Patient name if applicable				
Patient date of birth if applicable				
Patient hospital number if applicable				
Date & time received				
Received by (2 staff members)				
Checking step\data	Yes/No/N/A/ Data	Checker Initials		Date & time
Shipping documentation received: Shipping log Returns documents	Yes / No / N/A Yes / No / N/A			
Qualified Person (QP) Release Certificate / Certificate of Conformance/ Certificate of Analysis available (delete as applicable)	Yes / No			
Location of transport container checked on manufacturer's portal <b>if applicable</b> <sup>1</sup> (must be expected receiving centre/location)	Yes / No / N/A			
Stem Cell Lab (SCL) ATMP Batch Processing Record is complete and signed for release <b>if applicable</b> <sup>1</sup>	Yes / No / N/A			
Temperature log received/ downloaded from manufacturer or SCL and transport/storage temperature within limits <sup>2</sup>	Yes / No			
Any visible damage to shipper?	Yes / No			

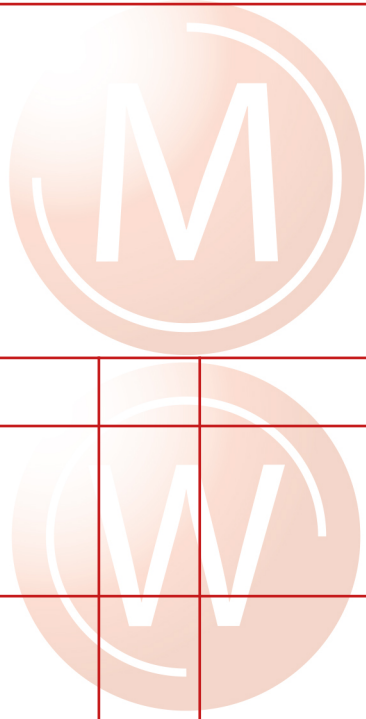Details and action taken				
Shipper ID matches delivery documentation	Yes / No			
Tamper-evident ties intact? Outer (Tag ID: _____) Inner (Tag ID: _____)	Yes / No			
Tamper-evident tie ID corresponds to documentation	Yes / No			
Transport container Data Logger within specification <sup>3</sup> (i.e. no alarms) on receipt?	Yes / No			
Quantity received (no. of bags/ vials/ syringes)	.....			
Number of bags/vials/syringes matches QP Certificate/ Certificate of Conformance/ Certificate of Analysis (circle appropriate parts)	Yes / No			
Product integrity visual check <sup>4</sup>	Pass / Fail			
Lot/Batch number				
Product batch/ID number matches QP certificate/ Certificate of Conformance/ Certificate of Analysis	Yes / No			
Name on product matches manufacturer portal / QP certificate / (circle applicable parts)	Yes / No			
Patient identifiers on product match manufacturer portal /QP certificate (circle applicable parts)	Yes / No			

Product dose matches Certificate of Conformance/ QP certificate/ Certificate of Analysis (circle applicable parts)	Yes / No			
Expiry date; Administration planned prior to expiry date?	..... Yes / No			
All documentation filed as per local policy, to be retained for period appropriate to product type?	Yes / No			
Comments:				

<sup>1</sup> ATMPs may be received by Pharmacy/clinical area direct from the manufacturer or via the local Stem Cell Laboratory. The checklist must be amended to reflect local processes. The batch processing records apply only to ATMPs which have been received by a SCL direct from the manufacturer for storage and then subsequently delivered to a clinical area.

<sup>2</sup> For some ATMPs, the temperature log for transport containers may only available in retrospect once the transport container is returned to the SCL/sending laboratory and data logger downloaded. Compliance with transit temperature parameters is assessed by the transport container data logger at the point of receipt.

<sup>3</sup> If temperature deviations or product defects have been identified, place the fresh cellular ATMP into appropriate storage at the correct temperature and label 'under quarantine'. Contact the manufacturer immediately and refer to local Standard Operating Procedures for ATMP product deviations. If an ATMP is leaking from its container, it should be appropriately stored as per local procedure to avoid contamination of other products or medicines and local spillage procedures must be followed.

<sup>4</sup> Product visual integrity check should be undertaken by checking for leakage and holding up the fresh cellular product to a light source. Further guidance on checking product integrity is detailed in the MW-ATTC training resource 'Visual inspection of cell therapy medicinal products', available via <https://www.theattnetwork.co.uk/resources>. If defects noted, refer to point 3.

## PLACING CELLULAR ATMP INTO STORAGE (if applicable)

Checking step\data	Yes/No/N/A Data	Checker Initials	Date & time
Required storage temperature range (from SPC or trial protocol or manual)			
Product placed into storage, storage location temperature verified and temperature constantly monitored?	Yes / No / N/A		
Refrigerator ID/Room ID (plus entered on local database if required)	Yes / No / N/A		
Receipt documented on manufacturers platform if required	Yes / No / N/A		

Completed receipt checklist sent to Pharmacy (if applicable)	Yes / No	Initials	Date & time

FINAL CHECK	Print name	Signature	Date
Checker 1			
Checker 2			